solac

MOD. ME7750

**MASAJEADOR PODAL
FOOT MASSAGER
MASSEUR DES PIEDS
FUßMASSAGEGERÄT
MASSAJADOR PODAL
MASSAGGIATORE PER I PIEDI
VOETMASSAGEAPPARAAT
MASÉR CHODIDEL
MASAŽER DO STÓP
MASÉR CHODIDIEL
LÁBMÁSSÁZS KÉSZÜLÉK
МАСАЖОР ЗА КРАКА
MASAŽER ZA STOPALA
APARAT PENTRU MASAJUL PICIOARELOR**



**INSTRUCCIONES DE USO
INSTRUCTIONS FOR USE
MODE D'EMPLOI
GEBRAUCHSANLEITUNG
INSTRUÇÕES DE USO
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING
NÁVOD K POUZITÍ
INSTRUCJA OBSŁUGI
NÁVOD NA POUZITIU
HASNÁLATI UTASÍTÁS
инструкция за употреба
UPUTE ZA UPOTREBU
INSTRUCȚIUNI DE UTILIZARE**

Solac is a registered Trade Mark

05/06

5- 7 • ESPAÑOL
8-10 • ENGLISH
11-13 • FRANÇAIS
14-16 • DEUTSCH
17-19 • PORTUGUÊS
20-22 • ITALIANO
23-25 • NEDERLANDS
26-28 • ČESKY
29-31 • POLSKA
32-34 • SLOVENSKY
35-38 • MAGYAR
39-41 • българск
42-44 • HRVATSKA
45-47 • ROMANA

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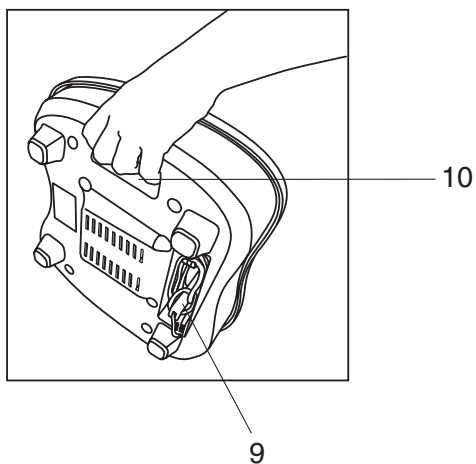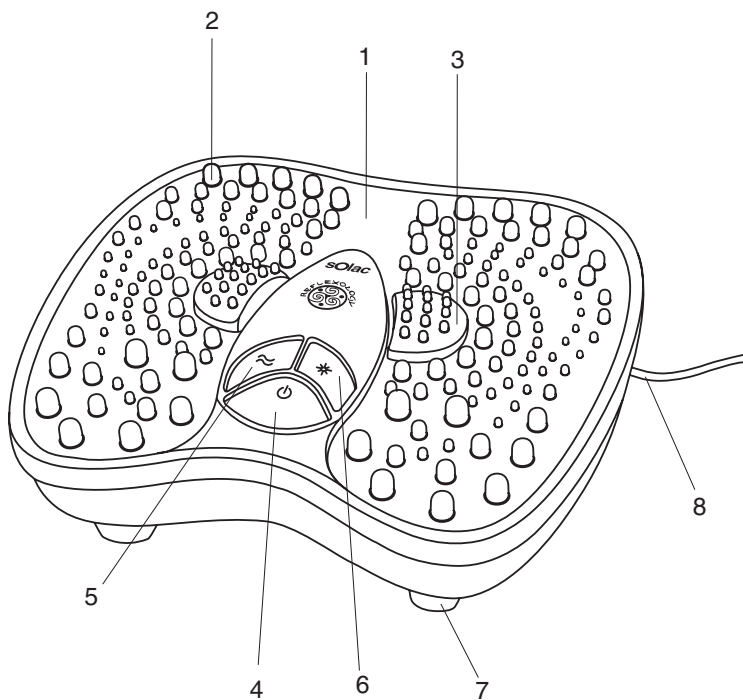FOOT MASSAGER

The reflexology massager by Solac has been developed in accordance with reflexology techniques to relax the body and improve general well-being. The digipuncture points act on the reflex areas of the body with a gentle, pleasant massage. The vibration of the massage points activates body circulation, releases tension and tones up foot muscles.

Thank you for purchasing our product. We hope you enjoy using your new Solac reflexology massager.

IMPORTANT

- Read the instructions carefully before using the appliance.
- Keep these instructions in a safe place for future reference.
- Before connecting the appliance, make sure that the voltage shown on the specifications plate is the same as the mains power supply at your home. Never connect it to a DC power supply.
- Remove all elements used for protection during transport or for sales promotions, such as paper or plastic bags, plastic sheets, cardboard and stickers, from the inside and outside of the appliance.
- Do not place the product in water or other liquids as it could be damaged.
- Keep the appliance out of the reach of children.
- Never unplug the massager by pulling on the power cable or by pulling on the product itself.
- Do not leave the power cable in contact with hot surfaces.
- Before using the massager, pull out the power cable fully and check that it is in perfect condition.
- Do not use the massager if the power cable is damaged, or the appliance is faulty. Take it to an Authorized Service Centre to have it repaired or replaced, as proper tools are required.
- Unplug the massager when it is not in use or for cleaning.
- Do not use accessories not recommended by Solac under any circumstances, as they could damage the appliance.
- This massager has been designed for home use. It can not be used on a continuous, uninterrupted basis.
- The appliance includes a heated surface. Persons sensitive to heat should take care when using the appliance.

- If you decide not to use the massager again, it is best to put it out of service. To do so, unplug the appliance and cut the power cable. In addition, it is best to disable the parts of the appliance that could be dangerous, especially for children, as they might use it improperly during play.

CAUTION! When you want to dispose of the appliance, NEVER throw it in the rubbish bin. Take it to your nearest CLEAN POINT or the waste collection centre closest to your home for processing. In this way, you will be helping to take care of the environment.

- Check the massager frequently for signs of wear or damage. If these signs appear or if the product has been used improperly, take it to an Authorised Service Centre.
- This massager has been designed for use when seated.
- Stop using the appliance if signs of irritation or discomfort appear.
- Electrodomésticos Solac, S.A. hereby declines responsibility for injuries or damage that could be caused to persons, animals or property due to the failure to heed these warnings.

MAIN COMPONENTS

1. Appliance body
2. Digipuncture points
3. Infrared heat points
4. Massager on/off button
5. Two-speed massage button
6. Infrared heat button
7. Non-slip legs
8. Power cable
9. Cable-winder
10. Carry handle

TIPS FOR USE

A) Duration of massage

The recommended duration of the massage is between 10 and 20 minutes. If you want a relaxing massage, the duration of the massage can be up to 20 minutes at minimum speed. If the massage is for therapeutic purposes, we recommend a maximum massage duration of 15 minutes, beginning with minimum intensity and changing to a higher-intensity speed after five minutes.

- b) If you suffer from any particular pathology or if you are pregnant, consult your doctor before using the reflexology massager.

- c) It is not recommendable to give very small children reflexology massages; rubbing the soles of their feet with your hand is sufficient. If you wish to give an older child a massage, use minimum speed and a duration of no longer than 10 minutes.
- d) We recommend using the reflexology massager with bare feet.
- e) It is very important to relax and sit comfortably for the massage. Sit comfortably and breathe deeply three times before beginning the massage.

INSTRUCTIONS

- Connect the appliance to the mains using the power cable (8).
- Sit comfortably and place your bare feet on the appliance. You will notice how your feet press gently on the digipuncture points that are going to give you the massage.
- Choose the type of massage you wish to receive:
 - If you prefer a gentle massage, press the massage button (5).
 - If you prefer a more intense massage, do not press button (5).
- To turn the massager on, press the on/off (4) button with your foot.
- If you wish to complete the massage with heat, press the heat button (6). The heat makes the massage more effective, activating the micro-circulation of the blood in the foot and toning up the muscles in the area.
- After use, disconnect the functions and turn off the appliance by pressing the on/off button again (4).
- Unplug the appliance from the mains by pulling gently on the plug.

MAINTENANCE

Clean the appliance with a damp cloth and, if necessary, use a very gentle detergent. Do not place the massager under water under any circumstances.

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GARANTE: ELECTRODOMESTICOS SOLAC S.A. Artapadura, 11 • 01013 VITORIA (SPAIN)