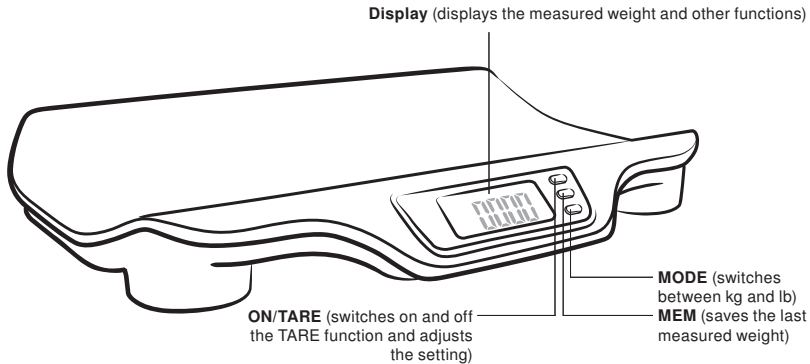
SENCOR®

SBS 4002

Before using this device, please read the user's manual thoroughly, even if you are already familiar with the use of other similar products. Use the device only as described in this user's manual. Keep this manual for future reference.

DESCRIPTION**SPECIFICATIONS**

Option to select between kg or lb
Maximum measurable weight: 20 kg/44 lb
Scale precision: 0.005 kg/0.01 lb
Memory for last measured weight
Low battery indicator (the display will show "Lo")
Error control (the display will show "Err")
Voltage: Four batteries type AAA
Automatic shut off

BUTTON DESCRIPTION

1. **ON/TARE**
Scale on
Function TARE (for example the weight of the child's towel)
Scale off. To switch the scale off, press and hold the button for about 3 seconds.
2. **MODE**
Changes the measuring units between kg/lb. By pressing the MODE button you switch between kg and lb.
3. **MEM**
Displays the last measured weight. If the scale shows 0,000, by pressing the button you may display the last saved value.

PREPARATION:

1. Battery Installation: Open the battery compartment cover on the bottom of the scale and place four AAA batteries inside. Observe the polarity.
2. Scale placement: Place the scale on even and flat surface.

NAKED CHILD WEIGHTING

1. Press the button ON/TARE to switch the scale on.
2. After few seconds the display will show 0.000 and the selected units.
3. If you wish to change the units, press the MODE button to select kg or lb. The selected units are shown on the right side of the LCD display.
4. Carefully place the child on the scale. After few seconds the scale will display the child's weight.
5. If you want to save the current weight (you can save the value flashing on the display), press the MEM button.

Note:

If the child is not steady during the weighting process, the correct weight may not be displayed. There will be no flashing value on the display. In such case, press the ON/TARE button shortly. The weight will be displayed while the button is depressed. WARNING. If you hold the button longer than 3 seconds the scale will switch off.

WEIGHTING A CHILD WITH TOWEL OR IN CAR SEAT ETC.

1. Press the button ON/TARE to switch the scale on.
2. After few seconds the display will show 0.000 and the selected units.
3. If you wish to change the units, press the MODE button to select kg or lb. The selected units are shown on the right side of the LCD display.
4. Place the towel or the car seat on the scale. After few seconds the scale will display the child's weight.
5. Press the ON/TARE button and the display will show 0,000.
6. Carefully place the child on the scale. After few seconds the scale will display the child's weight.
7. If you want to save the current weight (you can save the value flashing on the display), press the MEM button.

Note:

If the child is not steady during the weighting process the correct weight may not be displayed. There will be no flashing value on the display. In such case, press the ON/TARE button shortly. The weight will be displayed while the button is depressed. WARNING. If you hold the button longer than 3 seconds the scale will switch off.

CLEANING AND MAINTENANCE

Clean the scale with soft cloth.

Do not submerge the scale in water or any other liquid.

Do not use abrasive cleaning solutions.

Store the scale in dry, dust-free place and protect it from ultraviolet rays.

The storage temperature range shall be from 16°C to 40°C.

WARNING

- The scale must be placed on solid and even surface to ensure safe operations.
- We recommend changing the batteries for new ones as soon as the symbol "Lo" is displayed.
- Maintain the scale clean to ensure correct operation.
- Remove the battery when the unit is not being used for a long period of time.
- Make sure to place the child in the middle of the tray to ensure precise weighting process.
- If a load above 20,5kg is placed on the scale, an irreversible damage may occur. Such damages are not covered by the warranty policy.

Notes:

Due to continuous product improvement, the technical specifications contained herein may be changed without prior notification of the customer.

Using this scale is not a substitute for a paediatrician.

If you have any doubts, contact your medical doctor.

The importer of this product is not liable for wrong use or interpretation of measured values.

Disposal of used batteries

Do not dispose old batteries into regular household waste. Instead, use the appropriate waste collection and recycling centres. The improper disposal of this type of waste may fall under subject of national regulations for fines.

Liquidation of used electric and electronic equipment

This symbol on a product, its accessories or packaging signifies that this product must not be treated like household waste. Liquidate the product by submitting it to the collection site for recycling of electrical and electronic devices. Independent waste collection systems for the disposal of electric and electronic products exist in the European Union states and other European countries. Make sure to dispose off these products accordingly to help eliminate potential risks to the environment and human health as these products may be harmful and hazardous. Recycling waste materials helps maintain natural sources - for this reason, please do not mix your old electric and electronic facilities with your household waste. For further information on recycling, please contact an environmental protection worker from your local (municipal or district) office, workers of waste disposal sites or employees at the store, where you purchased this product.

Handling instructions for the disposal of used packaging material

Dispose off used packaging material at a designated site for waste disposal.

This product complies with the European Electromagnetic Compatibility Act.

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