solac

MOD.ME7710

Masajeador antcelulítico
Anti-cellulite massager
Masieur anticellulite
Anticellulite-Massagegerät
Massajador antcelulítico
Massaggiatore anticellulite
Massageapparaat tegen cellulitis

INSTRUCCIONES DE USO
INSTRUCTIONS FOR USE
MODE D'EMPLOI
GEBRAUCHSANLEITUNG
INSTRUÇÕES DE USO
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING

Solac is a registered Trade Mark

03/10

es	5-10 •	ESPAÑOL
en	11-16 •	ENGLISH
fr	17-22 •	FRANÇAIS
de	23-28 •	DEUTSCH
pt	29-34 •	PORTUGUÊS
it	35-40 •	ITALIANO
nl	41-47 •	NEDERLANDS

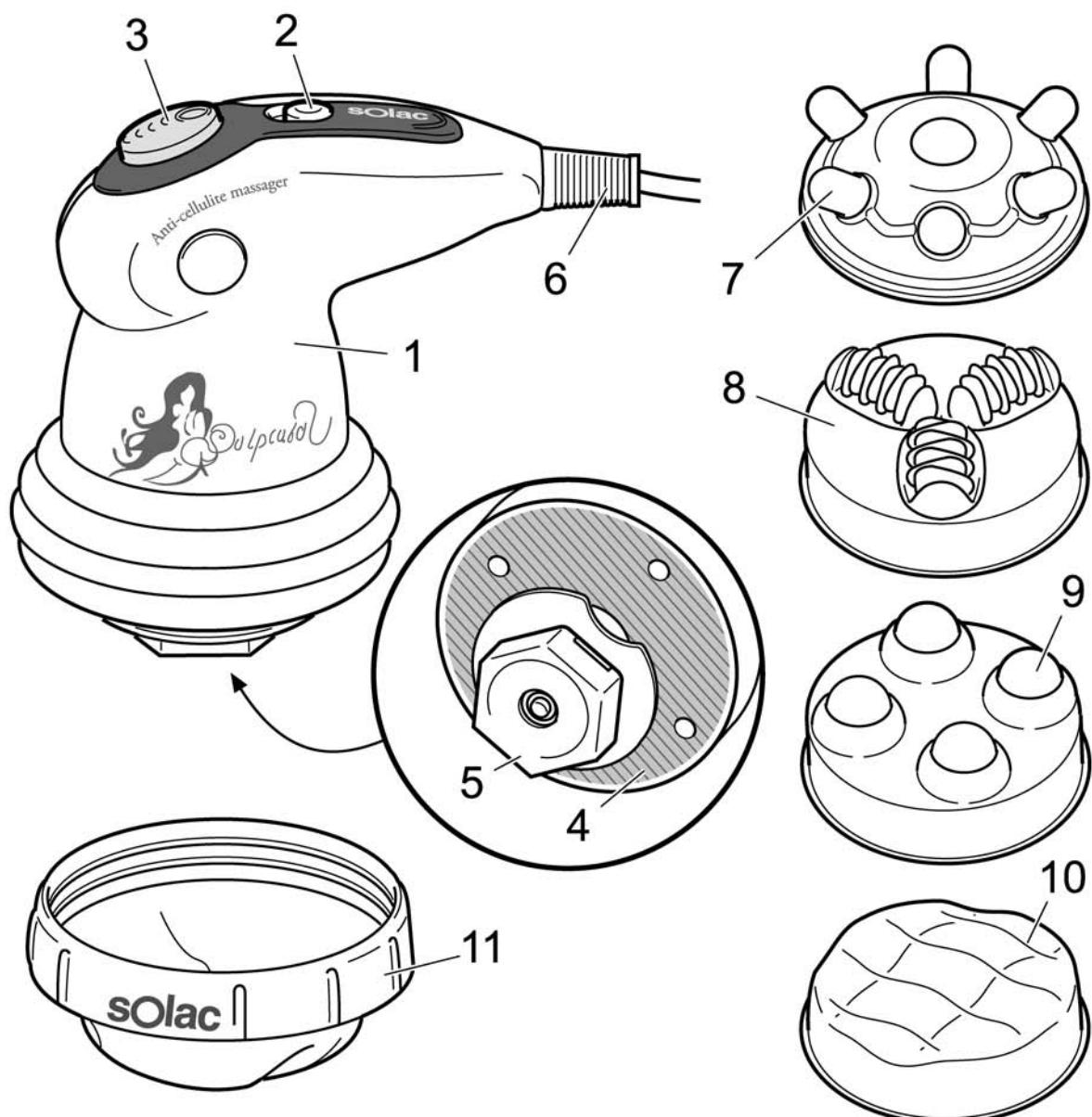
Fig. 1







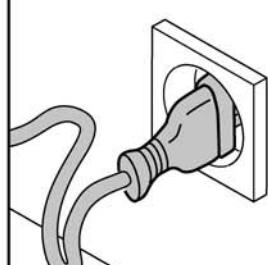
Fig. 2

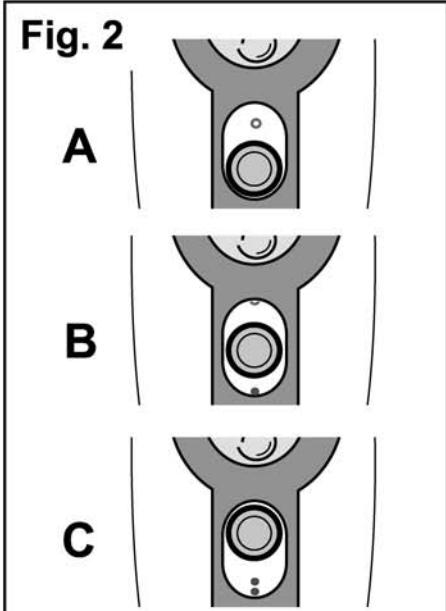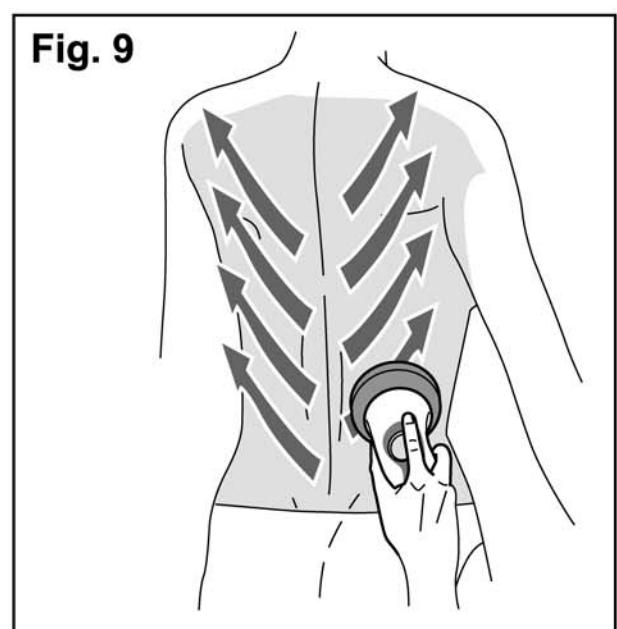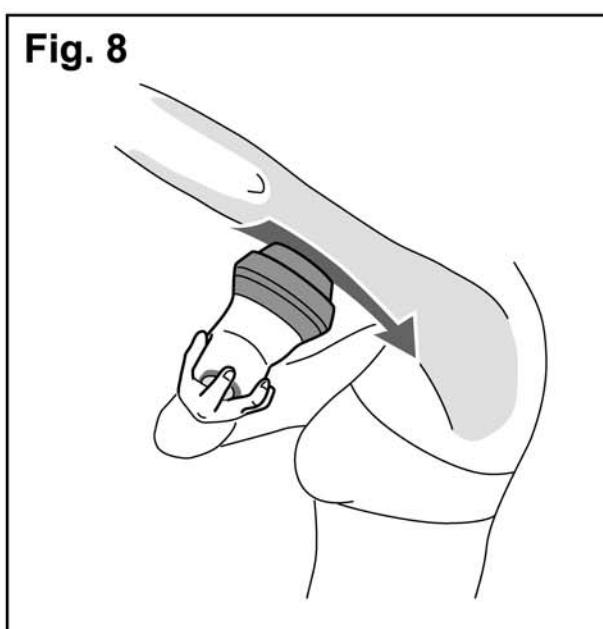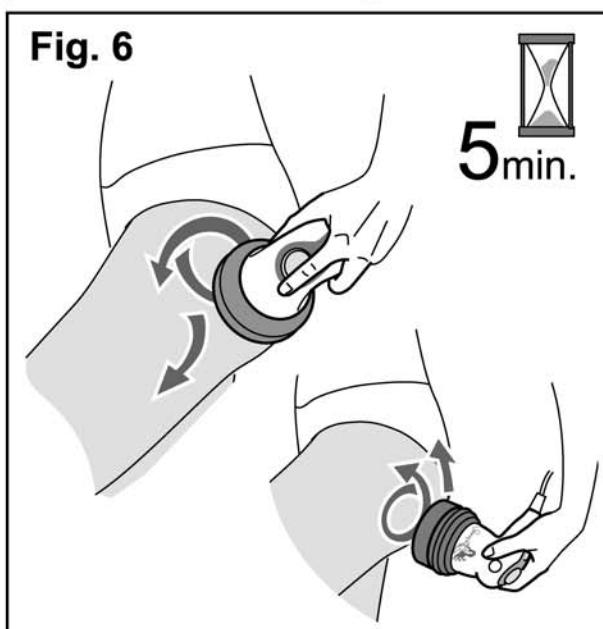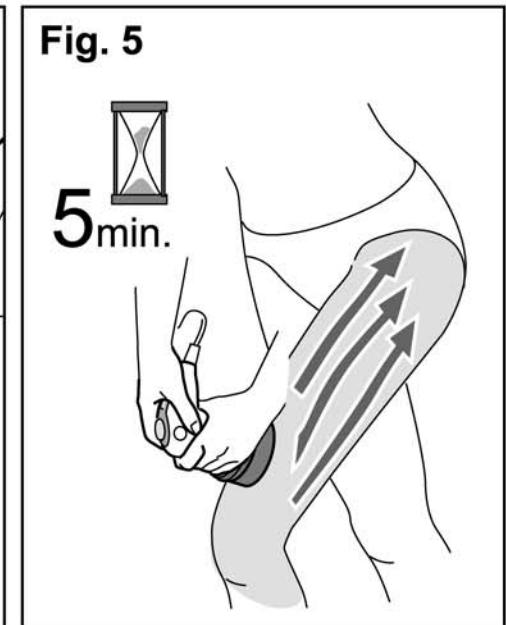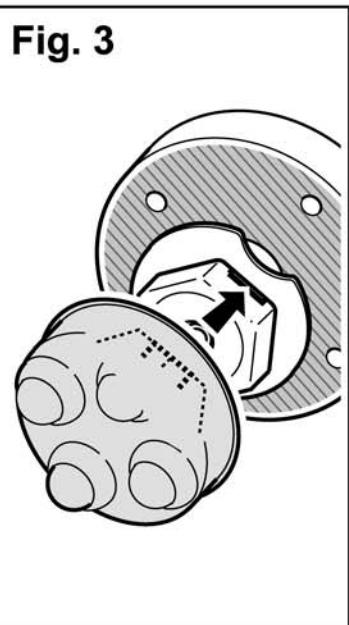. ESPAÑOL

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. ENGLISH

1 Safety Instructions

- Read these instructions completely before using your appliance.
- This appliance is only for domestic use, not professional use. Using it in any way other than according to the instructions could be hazardous. Do not use this appliance as a substitute for appropriate medical care.
- This appliance is not designed to be used by people (including children) whose physical, sensory or mental capacities are impaired or who lack experience or knowledge, unless they are supervised or have received instructions on how to use the appliance by a person who is responsible for their safety.
- Children must be supervised to ensure that they do not play with the appliance.
- Ask your doctor if you suffer from any health problem, for example, if you use a pace maker.
- Do not modify or repair the appliance. In there is any problem with the cord or any other part of the appliance, do not use it and take it to an authorised service centre.
- If the power cord is damaged, it must be replaced by the manufacturer or the manufacturer's after-sales service or by similarly qualified personnel to prevent possible hazards.
- Do not use parts or accessories not supplied or recommended by SOLAC.
- Before plugging in the appliance to the mains, check that the indicated voltage matches the voltage in your home.
- Always unplug the appliance after use and before removing or installing parts and before performing any maintenance or cleaning task. Also unplug it in the event of a power outage. When you plug in or unplug the appliance, it must be turned off.
- Do not pull on the power cord or hang the appliance by the cord. Always unplug the appliance by pulling on the plug, never the cord.

en

Do not use near bathtubs, sinks, showers or other recipients containing water. If it falls into water, do NOT try to pick it up. Unplug it immediately.

- Never leave the appliance exposed to the elements.
- Do not touch the appliance with wet hands when in operation.
- Never leave the appliance working on a wet surface or cloth.
- Do not use the massager if you are suffering from any skin rash or inflammation.
- **Do not use infrared radiation if you suffer from varicose veins, as it could exacerbate the circulatory problem.**
- **Do not use the massager in the abdomen area if you are pregnant or believe that you might be pregnant.**

IMPORTANT: When you want to dispose of the appliance, NEVER throw it in the rubbish bin. Take it to the CLEAN POINT or waste collection centre closest to your home for processing. You will thus be helping to care for the environment.

2 Description

Main parts

- 1 Main body of the appliance
- 2 Off/massage/massage + IR switch
- 3 Variable massage intensity control
- 4 Infrared radiation emitter
- 5 Massage head adaptor
- 6 Cord
- 7 Finger massage head
- 8 Roller massage head
- 9 Ball massage head
- 10 Wave-shaped head
- 11 "Wave-shape" head protector

Electrical characteristics

Charger power supply voltage: 230-240V, 50Hz, 22W

Professional Lipomodelling system.

Fights cellulite and flaccidity efficiently at home

This latest-generation massager allows you to enjoy an effective therapeutic massage at home and is especially recommended for people suffering from cellulite, fluid retention and flaccidity. Continuous use of the apparatus will help you to:

- Improve blood circulation.
- Obtain firmer and smoother skin.
- Eliminate fluids and toxins.
- Eliminate accumulated fat and "orange peel" skin.

3 Selecting the Massage Head

Head Types and Functions

- "Finger massage" head (7)
For painful cellulite and muscular pains.

- "Roller" head (8)
This eliminates cellulite, works against flaccidity, stimulates blood flow and return, tones the muscles, performs lymphatic drainage, reduces volume, relieves back ache, sciatica and lumbago, warms up the muscles, relaxes and relieves stress, for example.

- “Massage balls” head (9)
Stimulates and tones the muscles, helping to firm the skin. It improves circulation in the abdominal area, eliminating localised toxins and improving the intestinal tract, preventing constipation.
- “Wave-shape” head (10)
Tones and improves blood circulation. For use in the thigh and buttock areas.
For improved comfort, cover the wave-shape head (10) with the protective fabric (11).

Types of Cellulite

Identify the type of cellulite you suffer from, to help you choose the most effective massage head.

- *Oedematous Cellulite: Feels spongy, it does not hurt and moves when touched or when you change position.*
- *Painful Cellulite: Pain is felt when the area is pressed.*
- *Cystic Cellulite: Grainy and hard to the touch. Usually accompanied by stretch marks. Its temperature feels lower compared to the rest of the body.*

en

Treatment Type

ZONE	CONDITION	HEAD	INTENSITY
1 st STEP			
	Oedematous Cellulite Flaccidity Cystic Cellulite		
2 nd STEP			
LEGS AND BUTTOCKS			
	Painful Cellulite Tiredness		
	Fluid Retention		

ZONE	CONDITION	HEAD	INTENSITY
	Oedematous Cellulite Painful Cellulite Flaccidity		
ABDOMEN	Cystic Cellulite		
	Constipation		
BACK	Back Pain, Sciatica and Lumbago, Stress, etc.		
			
ARMS	Oedematous Cellu- lite Painful Cellulite Flaccidity		
	Cystic Cellulite		

4 General Recommendations for Use

- The intensity control allows you to regulate the intensity of the massage as required. Start at a gentle speed and progressively increase it, depending on your sensitivity to the massage.
- The infrared radiation is an additional function to the massage action. Use it to strengthen the effectiveness of the massage. This function moves fat, heating and draining it, as well as smoothing bumpy "orange peel skin".
- If using special treatment creams (firming and anti-cellulite creams, for example), apply these after the massage, which prepares the skin and

en

its tissues to absorb the ingredients.

You can use the "wave-shape" head (10) to apply the cream.

- Do not apply creams when using the protective fabric (11) on the "wave-shape" head.

5 Recommended Treatment Duration

- Apply a 2-3 minute massage on the abdomen and a 5 minute massage on the legs and buttocks. You may establish the duration of the massage depending on your sensitivity. The total treatment must be approximately 15 minutes.
- Do not use the appliance for more than 15 minutes. Using the appliance continuously for more than 15 minutes may cause it to overheat, reducing its life. Let it cool down before continuing use.
- Use Sculptural every day until you notice the results, which may be subsequently maintained by using it three days a week. To maintain the results over the long-term, adapt your lifestyle habits and food.

6 Use

☞ This appliance has been designed to provide therapeutic massages to help prevent cellulite, fluid retention and flaccidity. Do not use it for any other purpose.

- 1 Plug in the appliance (fig. 1).
- 2 Move the switch (2) to the massage position (fig.2 B) or to the "massage + infrared position (IR) position" (fig.2 C) according to the desired treatment.
- 3 Fit the head (7), (8), (9) or (10) corresponding to the desired treatment onto the head adaptor (5) (fig. 3).

☞ To choose which massage head to use, follow the instructions in the "Treatment Type" table.

- ☞ You may alternate the use of different heads during the same massage.**
- 4 Set the intensity of the massage using the intensity control (3) and place the massager on the area on which you are going to work.
 - 5 Set the intensity of the massage using the control (3), according to the area to be treated and its sensitivity to massage (see the "Treatment Type" table).

To do so, turn the control to one side or another with your finger, depending on the desired intensity (fig. 4).

- 6 Hold the massager firmly with both hands and move it over the skin in the direction of the blood circulation (see "Massage Direction").

- ☞ Without pressing too hard, simply glide the massager over the area to be treated.**
- 7 When you have finished using the massager, move the switch (2) to the "off" position (fig.2 A).

- 8 Unplug the appliance.

7 Massage Direction

Always glide the massager over the skin in the direction of the blood circulation.

- *Legs*

Place the appliance on the thigh, close to the knee, and glide it gently upwards towards the hip (fig. 5).

- *Buttocks*

Combine upwards and circular movements. You can end an upward massage on the legs with downwards and circular movements on the buttocks (fig. 6).

- *Abdomen*

Place the appliance on the navel and move it gently outwards (fig. 7).

 *Never make circular movements around the navel.*

- *Arms*

Place the appliance close to the elbow and glide it gently towards the armpit. Repeat the massage from the armpit towards the elbow (fig. 8).

- *Back*

Place the appliance in the area of the spine and glide it gently upwards towards the shoulder blade (fig. 9).

 *NEVER apply a massage directly on the spine.*

 *NEVER apply a massage on a bone (knee, ankle, elbow). The protruding parts of the massage head could knock the bone, causing an injury.*

8 Cleaning and upkeep

 *Do not use solvents, detergents or abrasive products.*

- 1 Unplug the appliance.
- 2 Use a dry cloth to clean the outside.

Storage

- When the appliance is not being used, unplug it and store it in a cool, dry place.
- To prevent the cable from becoming damaged, do not twist it, pull it or wind it around the appliance. If it is damaged, have it replaced immediately by qualified personnel.

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GARANTE: ELECTRODOMESTICOS SOLAC S.A. Artagan, 11 - 01013 VITORIA (SPAIN)